

# USER INSTRUCTIONS

## EASYLIFT HOME

### OVER-BED / OVER-CHAIR TABLE

S13880 Easylift Home Over-Bed Table

S13892 Easylift Home Over-Bed / Over-Chair Table

#### ⚠ General Safety Warning ⚠

Please make sure that you read and follow these instructions carefully. Keep these instructions in a safe place for future reference.

## 1. INTENDED USE

To provide users with a firm, flat surface that is suitable for eating, writing, reading, etc. The single column construction and adjustable height makes them suitable for use with most types of beds and chairs. This product is NOT suitable for riser recliner chairs.



S13880



S13892

### S13880: Easylift Home Over-Bed Table

Suitable for most divan and slatted beds, as well as height adjustable medical beds.

### S13892: Easylift Home Over-Bed / Over-Chair Table

Suitable for most divan and slatted beds, as well as height adjustable medical beds. It is also suitable for use with most chairs (excluding riser recliner chairs).

## 2. SAFETY



**WARNING:** The central column is spring loaded. Make sure the central column hand-wheel is in position & fully tightened before assembly/ disassembly of the table.

**WARNING:** Do not exceed the safe working load of 25 kg (4 st).

**ALWAYS** make sure that the table is correctly positioned and adjusted to allow it to be used safely and comfortably.

**DO NOT** lean or sit on the table.

**DO NOT** place hot objects directly on the table surface. Maximum heat resistance 65°C (149°F).

**MAKE SURE** the area is clear before moving the table.

**TAKE CARE** when moving the table over rough or uneven surfaces. Lift the table over large ridges and wide gaps.

**ENTRAPMENT HAZARD:** Do not use this product within 318 mm of headboard or other similar object.

**ENTRAPMENT HAZARD:** Always move the table to a safe location when the user is sleeping.


## Safe Working Load



## 3. UNPACKING

Remove the packaging and check all parts carefully. Contact your supplier immediately if parts are missing or damaged and do not try and assemble.

### 3.1 Contents:

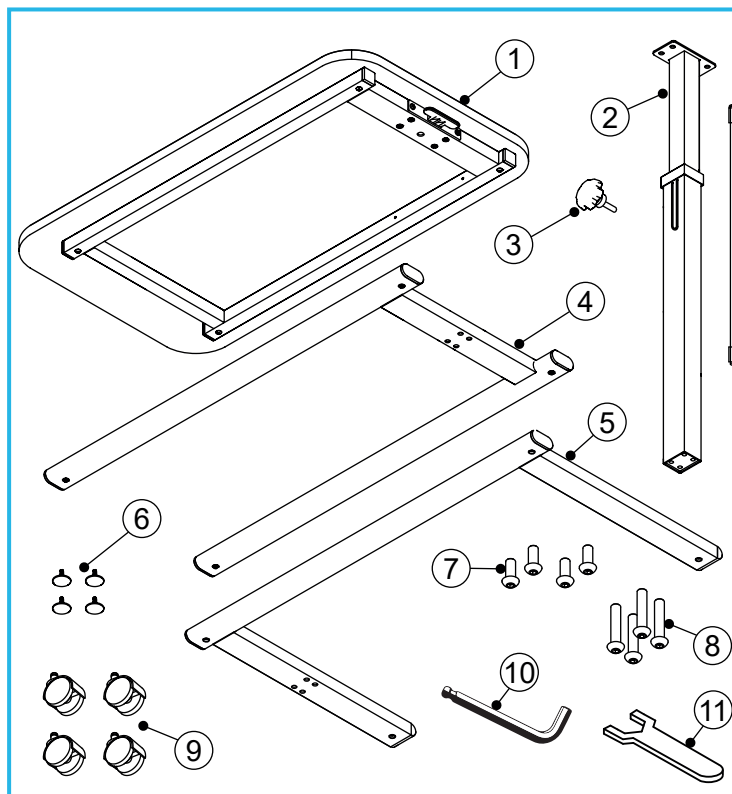


1. Table Top Assembly x 1
2. Central Column Assembly & Inner Rod x 1  
**Note:** Do not remove the PVC anti-rattle sleeves from either end of the inner rod.
3. Central Column Handwheel x 1
4. Base Assembly x 1 (S13880 only)
5. Base Assembly x 1 (S13892 only)
6. Mushroom Feet x 4
7. M8 x 20 mm Socket Head Button Screws x 4
8. M8 x 40 mm Socket Head Button Screws x 4
9. Castors x 4
10. Allen Key x 1
11. Castor Spanner x 1

## 4. ASSEMBLY


**Note:** We recommend a 2 person assembly.

### 4.1 Fit the Central Column Handwheel

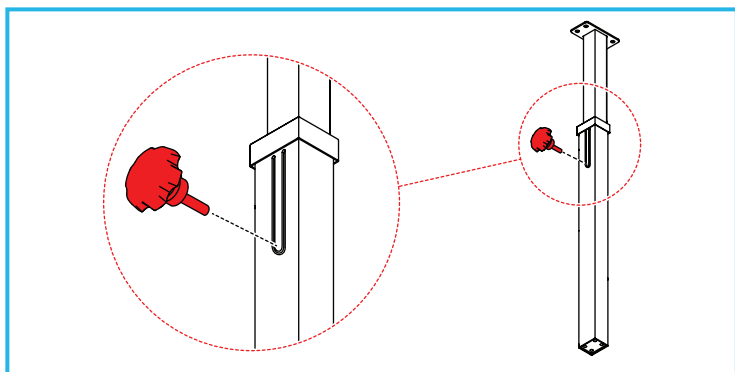


Insert the hand-wheel and turn clockwise to tighten.

### 4.2 Assemble the Table Top & Central Column (S13880 & S13892)

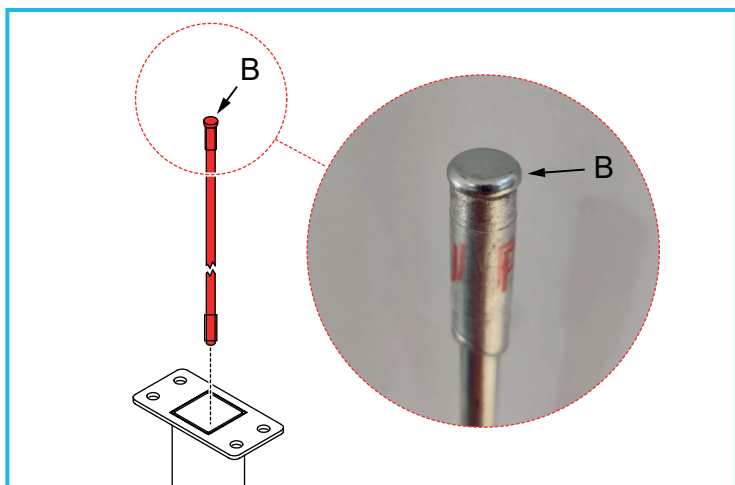

1. Slide the inner rod (B) into the central column.

**Note:** Make sure the flanged end is facing the correct way up. The plastic anti-rattle cover is marked 'UP'.

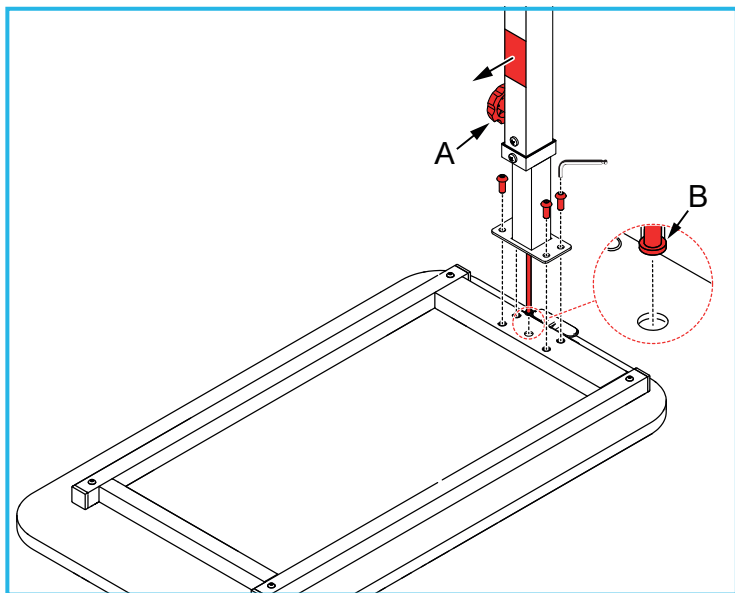

2. Guide the flanged end of the inner rod (B) into the height adjustment lever hole in the main frame.


**Note:** Make sure the label on the central column is facing inwards before assembly.

3. Secure the column assembly to the table top assembly using the 4 x M8 x 20 mm socket head screws and the Allen Key provided.

4. Make sure that all four screws are fully tight before you continue the assembly.

### 4.3 Assemble the Base to the Central Column (S13880 only)


**Note:** Make sure that the base faces the same direction as the table top and the castor inserts face upwards (as shown).

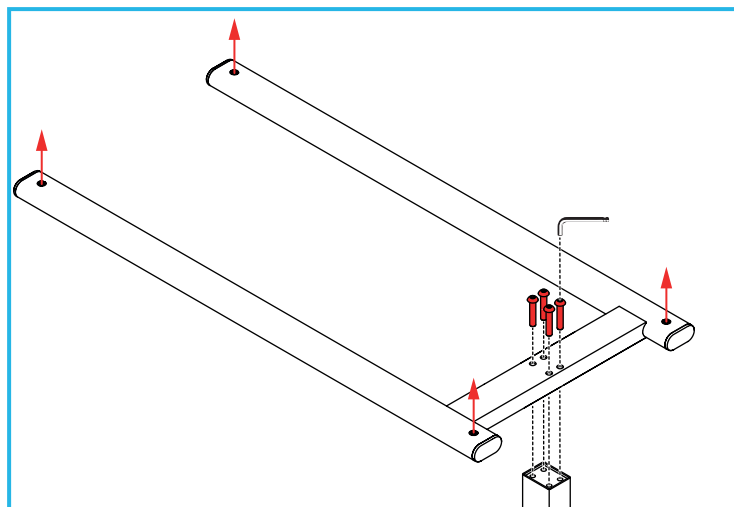

1. Secure the base assembly to the column assembly using the 4 x M8 x 40 mm socket head screws and the Allen Key provided.
2. Make sure that all four screws are fully tight before you continue the assembly.

### 4.4 Assemble the Base to the Central Column (S13892 only)

**Note:** Make sure that the base faces the same direction as the table top and the castor inserts face upwards (as shown).

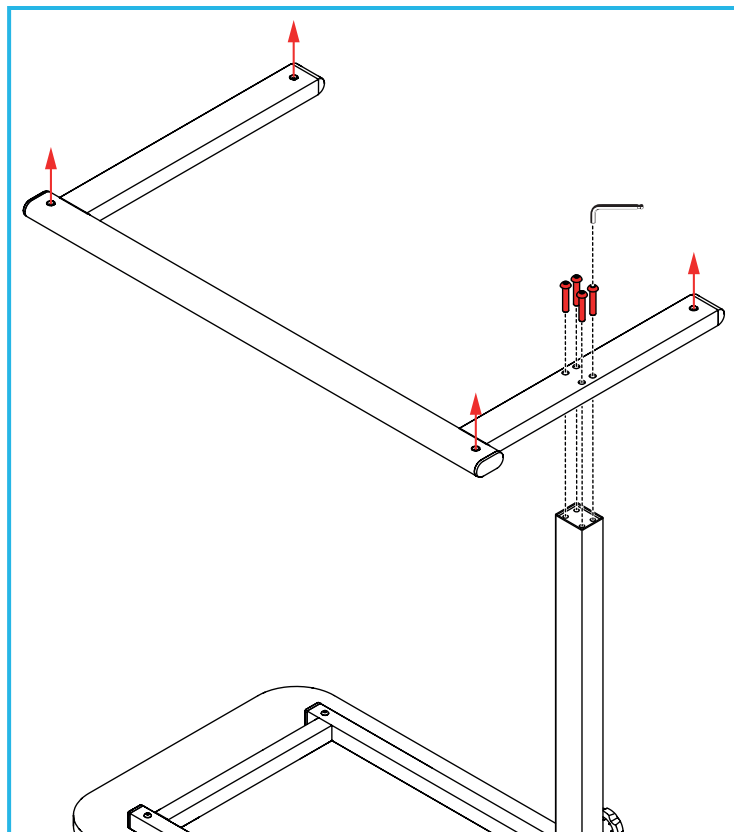


1. Secure the base assembly to the column assembly using the 4 x M8 x 40mm socket head screws and the Allen Key provided.

2. Make sure that all four screws are fully tight before you continue the assembly.

## 4.5 Fitting the Castors / Feet

### 4.5.1 Fitting the Castors (S13880 & S13892)

**Note:** The castor fitting method is the same for both products, S13880 assembly is shown in the figure below.

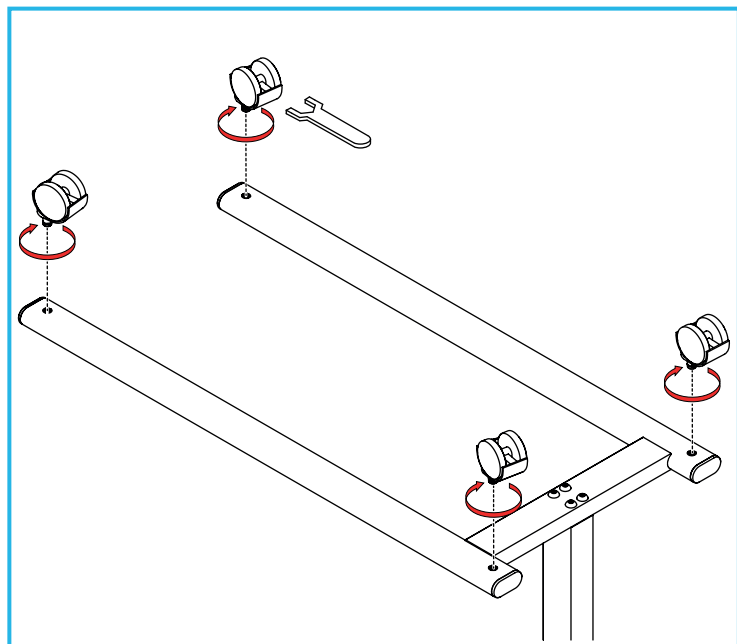


1. Insert the screw thread into the threaded insert, hold the castor and turn the nut by hand until the threads engage.
2. Tighten the castor threads using the 14mm spanner provided.
3. To remove the castors reverse step 1 & 2.

### 4.5.2 Fitting the mushroom feet (S13880 & S13892)

The mushroom feet should be fitted for use with a low base divan.

**Note:** The mushroom feet fitting method is the same for both products, S13880 is shown in the figure below.

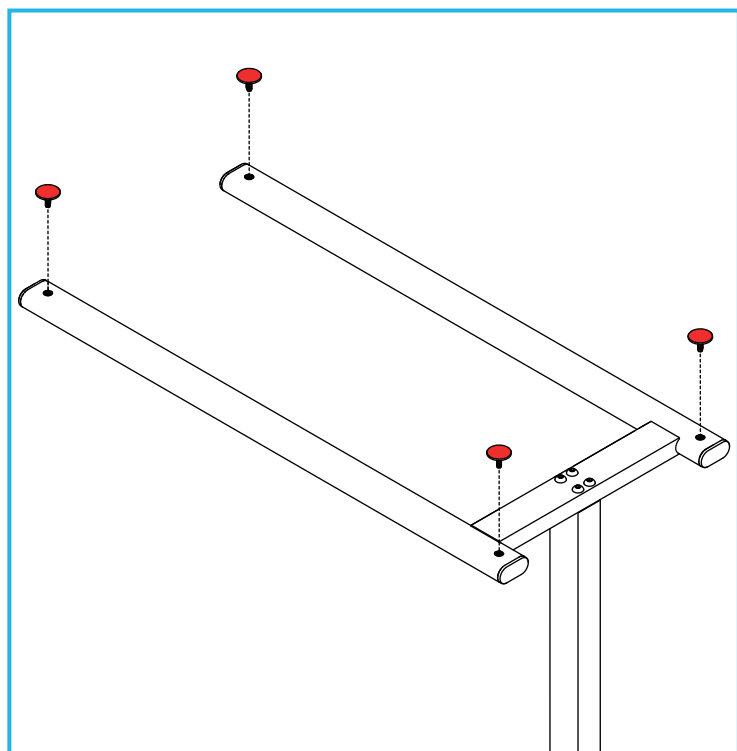


Press-fit the mushroom feet into the castor insert holes.


## 5. IN USE

### 5.1 Locking and Unlocking the Brakes.

The table has 4 independently lockable castors. We recommend that the castor locks are always fully engaged when the product is in use or stationary.

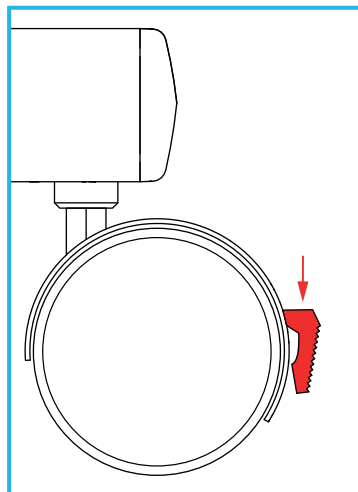


**Fig 1: Castor Unlocked**

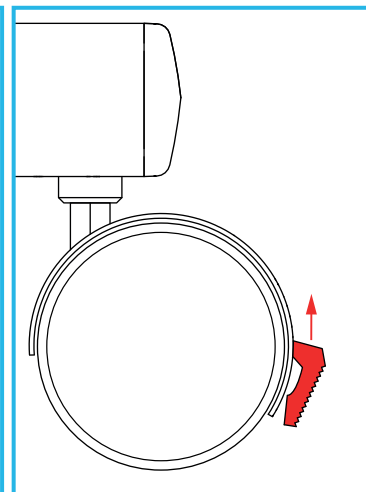

**Fig 2: Castor Locked**

1. To lock the castor brake, press down on the castor tab (Fig 1).
2. To disengage the castor brake, flip up the castor tab (Fig 2).

### 5.2 Raising & Lowering the Table

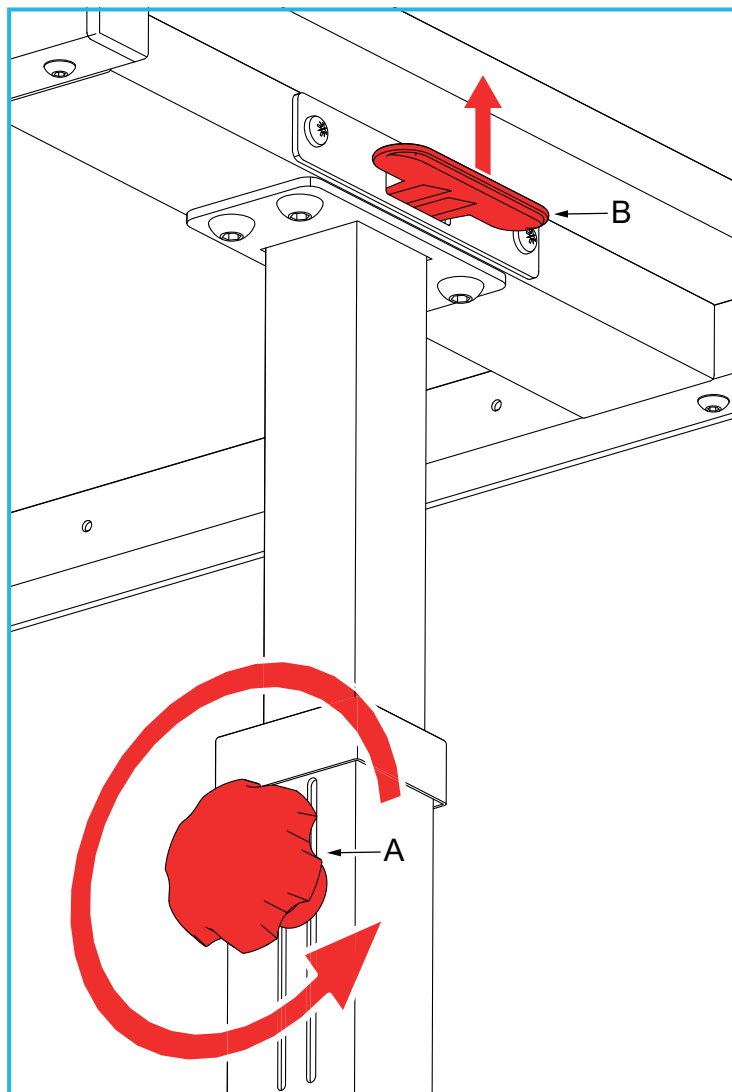

1. To release the central column, unscrew the locking handwheel (A) by turning anti-clockwise.

2. Pull lever (B) up to release the spring mechanism, then simultaneously raise or lower the table-top into position.

3. To lock the central column, screw the locking hand-wheel (A) clockwise.

## 6. MAINTENANCE

Regularly check the wheels rotate freely.

Regularly check the brakes function correctly in accordance with Section 5.1.

Always check the locking handwheel is fully tight before use.

Make sure the hand-lever is functioning correctly.

## 7. CLEANING INSTRUCTIONS

Make sure that all dirt and debris is removed from around the wheels.

Clean regularly using a general purpose, neutral detergent and warm water.

Please check the suitability of the cleaning products before use. Do not use abrasive cleaners or cloths as these could permanently damage the surface and invalidate any warranty.

DO NOT use solvents to clean the product (if painted)

To disinfect, wipe with a dilute bleach solution, and rinse thoroughly with fresh water.

Allow product to fully dry prior to use/storage.

## 8. ACCESSORIES & SPARE PARTS

S20161: Spares Kit

- Replacement Mushroom Feet (set of 4).
- Replacement Castors (set of 4).

For more details please visit our website or call Customer Services on the details below.

## 9. PRODUCT DISPOSAL

This product must be thoroughly cleaned/sterilised prior to disposal.

This product is manufactured from metal, wood & plastic and can be disposed of through Local Authority Household Waste and Recycling Centres or through Local Trade Waste Collectors.

## 10. PRODUCT GUARANTEE

NRS Healthcare guarantee this product for a period of 24 months, from the date of purchase, against faulty materials and workmanship.

This guarantee does not cover abuse, interference or tampering with any part of the products. Issues arising from fair wear & tear and customer misuse are not covered.

This guarantee covers domestic use only. If the product is used in commercial or semi-commercial settings such as care homes, schools, rental or communal accommodation etc, the guarantee is valid for one year only. Consumable parts are also guaranteed for one year only.

Any guarantee issued by NRS Healthcare does not affect your statutory rights.

**NRS HEALTHCARE** Sherwood House, Cartwright Way, Forest Business Park  
Bardon Hill, Coalville, Leicestershire, LE67 1UB, UK.

Tel: Customer Services: +44 (0)345 121 8111

Email: [customerservice@nrshealthcare.co.uk](mailto:customerservice@nrshealthcare.co.uk)

NRS Healthcare is a leading supplier of community equipment and daily living aids.  
Please visit one of our websites to view our full range of products:  
[www.healthcarepro.co.uk](http://www.healthcarepro.co.uk) or [www.completecareshop.co.uk](http://www.completecareshop.co.uk)