

User Instructions

Overbed/Over Chair Tables



Overbed/Over Chair Table
- Beech and White Frame S37779



Overbed/Over Chair Table
- Dark Oak and Black Frame S37781











Read and follow these instructions carefully. Failure to do so can result in injury. Keep these instructions in a safe place for future reference.

Intended Use



Used for working, eating or reading while in a chair or bed. The dual column design means that it can work as a flat surface or tilt in either direction, and works on either side of the bed or chair.

Symbols


Symbols used on the product and in these instructions.

	Manufacturer		Batch Reference Code		Consult the instructions for use
	Item/Catalogue number		Important Note		Recycling Symbol
	Maximum Load Weight		Warning		

Safety



Do not exceed the maximum load weight of 10 kg (1.6 st).



Inspect the product before use. Contact the service provider if you see any damage.

MAKE SURE

- The table is correctly positioned and adjusted to allow it to be used safely and comfortably.
- The area is clear before moving the table.

TAKE CARE when moving the table over rough or uneven surfaces. Lift the table over large ridges and wide gaps.

DO NOT:


- Use in wet areas such as bathrooms or wet-rooms.
- Stand or sit on the table.
- Use the table as a walking aid or to transport items.
- Put hot objects directly on the table surface. Maximum heat resistance 65°C (149°F).

ENTRAPMENT HAZARD:

- Always move the table to a safe location when the user is sleeping.

Unpacking

Remove the packaging and check all parts carefully.






Contact your supplier immediately if parts are missing or damaged. Do not try to assemble.

Contents

	<div>1. Table top</div> <div>2. Telescopic Tubes x 2</div> <div>3. Castors x 4 (2 with brakes and 2 without brakes)</div> <div>4. E-Clips x 2</div> <div>5. Vertical tubes x 2</div> <div>6. Slider feet (for use with divan bed)</div> <div>7. Base Frame</div> <div>8. Spanner</div>
	<div>3. Castors with brake x 2</div>
	<div>3. Castors without brake x 2</div>

Assembly/Installation

1

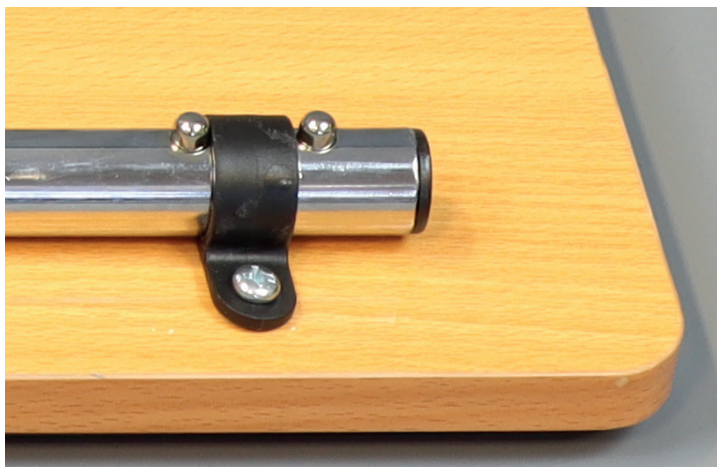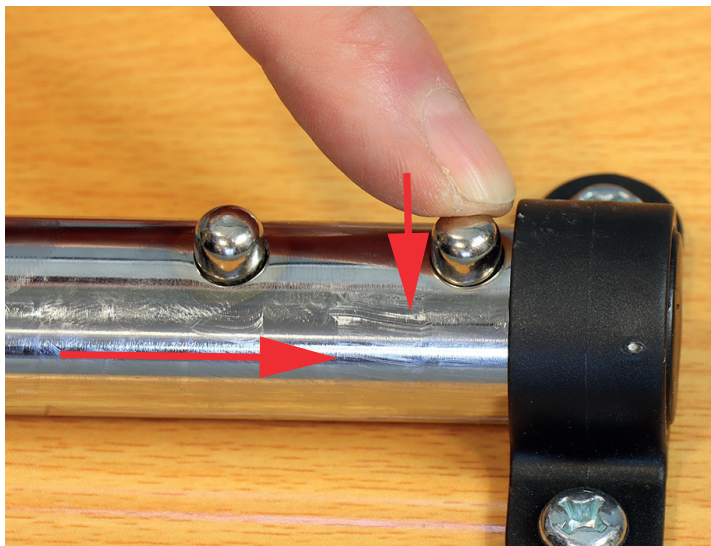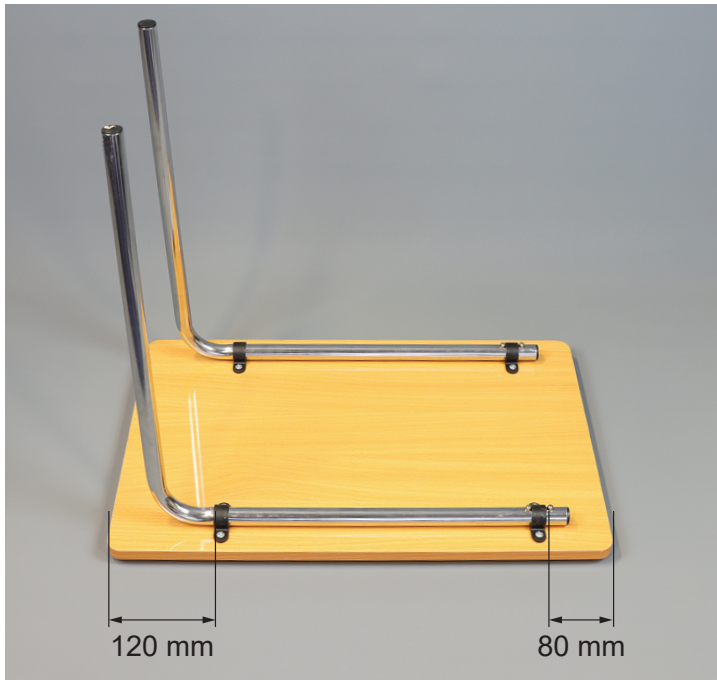

Put the table top upside down on a clean surface.

Insert the Telescopic Tubes into the saddles, pressing the button clips to allow the tubes to pass through the saddles.




The fittings under the table top are not equal distances from the edge. See image.

2



Decide whether you need to fit castors or the sliders for use with a divan bed.

Use the spanner to attach the four castors to the threaded holes in the Base Frame. Fit the two castors with brakes to the holes near the ends of the frame.

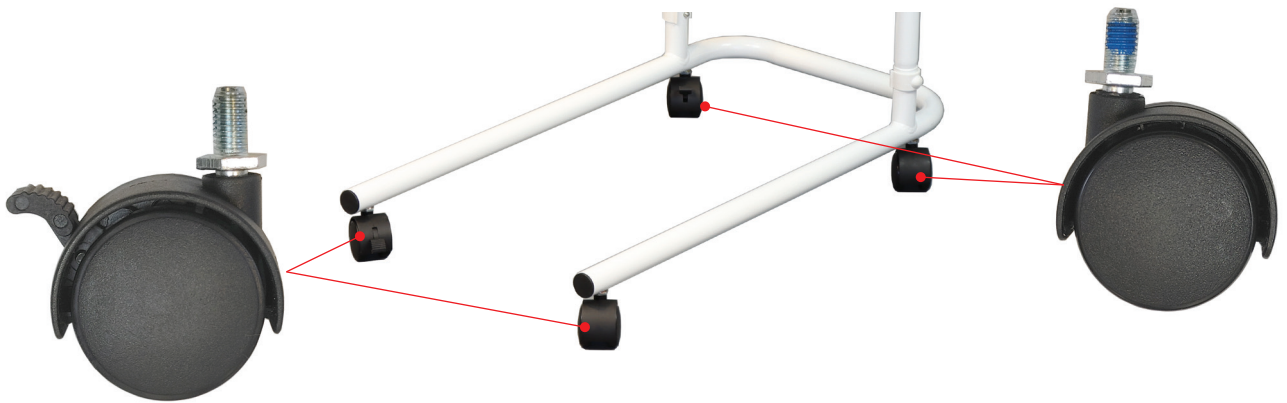

Or, fit the four sliders by simply pressing them into the threaded holes.

3



Use the E-clips to fasten the vertical tubes to the Base Frame.



Fix the vertical tubes so that the handwheels face the outside of the product.

4



Undo the handwheels and insert the telescopic tubes into the vertical tubes.


5





Table assembled.

Adjusting and Using the Table

1

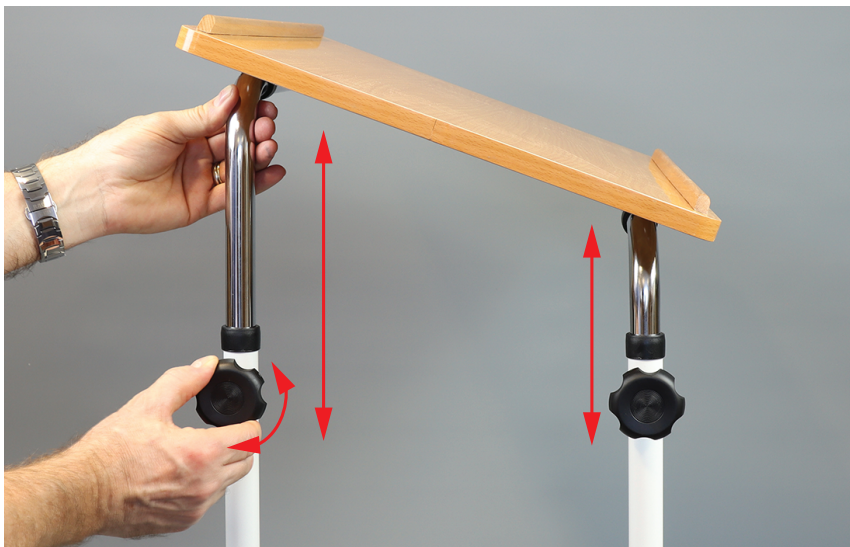


Undo and re-tighten the handwheels one at a time to adjust the height and angle of the table.

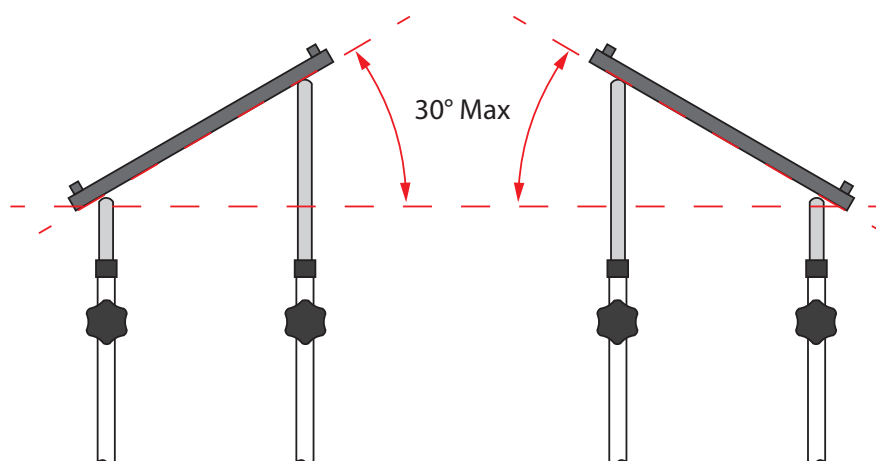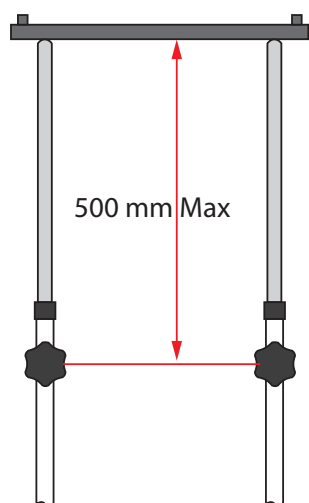

Adjust the table within the maximum height and angle limits.

2

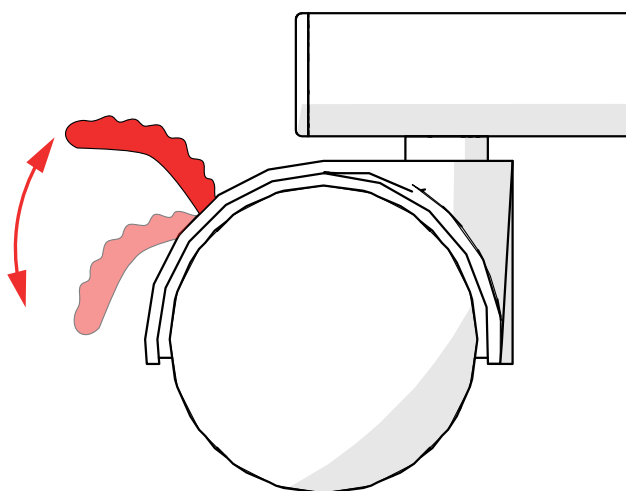

Lift up the small tab to unlock the castors if fitted.
Press down the tab to lock the castors.

Maintenance

Inspect the product at regular intervals and before re-use for signs of general wear or damage.
Report to the supplier any suspected faults or damage and take the product out of use immediately.
Check the fixings at regular intervals.

Cleaning

Clean regularly using a general purpose, neutral detergent and warm water. Check the suitability of the cleaning product before use. Do not use abrasive cleaners or cloths as these could permanently damage the surface and invalidate any warranty.

To disinfect, wipe with a dilute bleach solution, and clean with a damp cloth.

Allow product to fully dry before using and storage.

Accessories & Spare Parts

- S28249 Spares (2 x E-Clips, 4 x Castors, 4 x Sliders)

For more information, see the end of this document for our contact details.

Product Disposal

Thoroughly clean this product before disposal and recycle where possible.

This product is made of wood, metal and plastic. It can be disposed of through Local Authority Household Waste and Recycling Centres or through Local Trade Waste Collectors.



Recycle any cardboard packaging.

Product Guarantee

NRS Healthcare guarantee this product for a period of 24 months, from the date of purchase, against faulty materials and workmanship.

This guarantee does not cover abuse, interference or tampering with any part of the products. Issues arising from fair wear & tear and customer misuse are not covered.

This guarantee covers domestic use only. If the product is used in commercial or semi-commercial settings such as care homes, schools, rental or communal accommodation etc, the guarantee is valid for one year only. Consumable parts are also guaranteed for one year only.

Any guarantee issued by NRS Healthcare does not affect your statutory rights.

If you have any problems with our products, please contact us. See the end of this document for our contact details.

We are a leading supplier of community equipment and daily living aids.
Please visit one of our websites to see our full range of products.

NRS Healthcare. Sherwood House,
Cartwright Way, Forest Business Park,
Bardon Hill, Coalville,
Leicestershire, LE67 1UB, UK

Tel: +44 (0)345 121 8111
Email: customerservice@nrshealthcare.co.uk
Web: www.healthcarepro.co.uk
or www.completecareshop.co.uk