

USER INSTRUCTIONS

ECONOMY FOLDING SHOWER STOOL

L97729

Economy Folding Shower Stool

⚠ General Safety Warning ⚠

Please make sure that you read and follow these instructions carefully. Failure to do so can result in injury. Keep these instructions in a safe place for future reference.

Maximum User Weight



1. INTENDED USE

For use in a shower to provide a stable, seated support, where a user is unable to stand for long periods (i.e. during showering). It folds easily for storage and is ideal for travel use. The textured seat with drainage holes helps avoid the possibility of slipping and is ideal for use when space is limited.



2. SAFETY

Warning: The maximum user weight must not be more than 136 kg (21 st).

Avoid sitting down heavily onto the seat.
Do not place on an uneven or sloping surface.
Make sure the seat is stable and flat before sitting.
Check product regularly for signs of wear or damage.
Do not use in baths.

3. UNPACKING

Remove the packaging and check all parts carefully. Contact your supplier immediately if parts are missing or damaged and do not try and assemble.

Contents:


Seat Assembly x 1

Leg Extensions x 4

ASSEMBLY/INSTALLATION



The leg extensions are numbered and must be inserted into the corresponding seat frame tube.

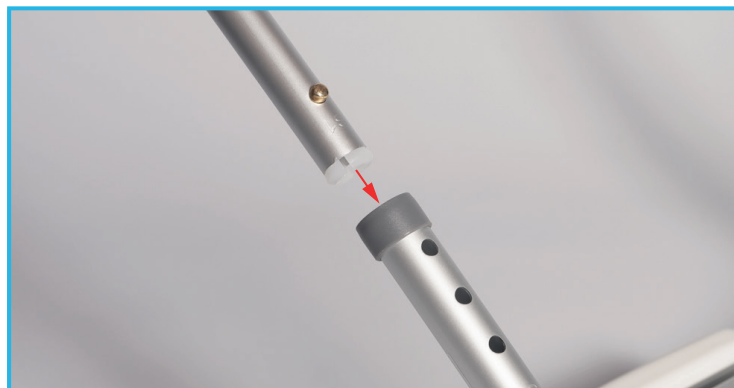

1. Inset the leg extensions numbered S1 into the seat frame tubes with the corresponding numbers (S1).
2. Inset Leg extensions numbered S2 into the seat frame tubes with the corresponding numbers (S2).



Make sure that the button clips have fully clicked into the adjustment hole and that all legs have been adjusted to the same height.



3. Push the seat frame support tube into the seat frame stop.



Make sure that the seat frame support tube has fully snapped into the seat frame stop.

4. IN USE

4.1 Using the Shower Stool

1. Place the assembled seat in the shower.
2. Make sure the seat is fully open and that the frame tube is securely located.
3. Make sure all four legs are level and in contact with the shower floor.



4. Before sitting down, check the seat is fully secure, by pushing down on the seat firmly with both hands.
5. Sit down centrally on the shower stool.

4.2 Folding the Shower Stool



1. Use your thumb and fingers to gently lever the seat frame out of the seat frame stop.



If the seat frame does not lift out of the seat frame stop, use your other hand in the same way to apply a little more pressure.

2. Fold the shower stool and store in a safe place.

5. MAINTENANCE

Inspect the stability and rigidity of the product on a regular basis.

Inspect the product regularly for wear and damage. Do not use if the rubber feet are damaged.

6. CLEANING INSTRUCTIONS

Clean regularly using a general purpose, neutral detergent and warm water.

Please check the suitability of the cleaning products before use. Do not use abrasive cleaners or cloths as these could permanently damage the surface and invalidate any warranty.

To disinfect, wipe with a dilute bleach solution, and rinse thoroughly with fresh water.

Allow product to fully dry prior to use/storage.

7. PRODUCT DISPOSAL

This product must be thoroughly cleaned/sterilised prior to disposal.

This product is manufactured from metal and plastic and can be disposed of through Local Authority Household Waste and Recycling Centres or through Local Trade Waste Collectors.

8. PRODUCT GUARANTEE

NRS Healthcare guarantee this product for a period of 24 months, from the date of purchase, against faulty materials and workmanship.

This guarantee does not cover abuse, interference or tampering with any part of the products. Issues arising from fair wear & tear and customer misuse are not covered.

This guarantee covers domestic use only. If the product is used in commercial or semi-commercial settings such as care homes, schools, rental or communal accommodation etc, the guarantee is valid for one year only. Consumable parts are also guaranteed for one year only.

Any guarantee issued by NRS Healthcare does not affect your statutory rights.

NRS HEALTHCARE Sherwood House, Cartwright Way, Forest Business Park
Bardon Hill, Coalville, Leicestershire, LE67 1UB, UK.

Tel: Customer Services: +44 (0)345 121 8111

Email: customerservice@nrshealthcare.co.uk

NRS Healthcare is a leading supplier of community equipment and daily living aids.

Please visit one of our websites to view our full range of products:

www.healthcarepro.co.uk or www.completecareshop.co.uk