



Care & Independence

Sling Model:

Serial Number:

User Name:

*The sling model and serial number
are located on the label on the reverse of the sling.*

Care & Independence Ltd is a UK manufacturer registered with MHRA as a Class 1 Medical Device Manufacturer and this sling is a Class 1 medical device as defined in Article 51, Annex VIII Rule 1 of the Medical Device Regulation 2017/745.

COMFORT VALUE™ slings are designed and manufactured to meet or exceed BS EN 10535:2006 for the transfer of disabled persons. COMFORT VALUE™ slings are designed and manufactured by Care & Independence, a certified ISO9001:2015 organisation.

Declaration

Care & Independence Ltd, declares under its sole responsibility that the aforementioned products to which this declaration relates, and which bear the CE marking, are in conformity with the applicable requirements of EU regulation 2017/745 as transposed in the laws of the Member States.

Intended Purpose

The sling is intended to be used with hoisting equipment to allow a care giver, or the individual needing assistance, to reposition a person, to a different standing, or sitting position, to allow numerous activities to be completed.

THIS USER GUIDE IS APPLICABLE
FOR THE FOLLOWING
COMFORT VALUE™ RECLINE SLINGS:

CVR: COMFORT VALUE™ Recline

CVR-HB: COMFORT VALUE™ Recline Highback

ISSUE DATE: 31/01/23
VERSION: UG/CVRECLINE-EN06
STATUS: APPROVED
AUTHORISED BY: Technical Director

.....
J. ASH

1. GENERAL



1.1. WARRANTY: The COMFORT VALUE™ slings covered in this User Guide have been awarded a THREE YEAR WARRANTY. Any manufacturing defect will be rectified free of charge. Product damage through wear and tear, neglect or accidental damage is not covered by the warranty.

1.2. COMPATIBILITY: Care & Independence slings have been manufactured to work with most of the major hoist manufacturers' equipment. A compatibility statement is available on request.

2. SLING SAFETY



2.1 GENERAL:

2.1.1. Only competent persons should ever use hoists and slings.

2.1.2 A risk assessment must be completed before using the sling, even if the hoist and sling is provided from the same manufacturer.

2.1.3 It is essential that the sling is fit for purpose; the correct size and type of sling for an individual is essential to ensure the safety and comfort of the person being transferred.

2.1.4 Product care and inspection instructions should be followed at all times.

2.1.5 **No sling must ever be used to lift more than its rated load** (see label).

2.1.6 The intended sling and hoist combination should be compatible.

2.1.7 **If there is any doubt about sling safety, remove the sling from use until checked and authorised for use by a competent person.**

2.2 SLING INSPECTION:

2.2.1 Care & Independence recommend pre-use sling checks by competent person.

2.2.2 Slings are Class One Medical Devices and are subject to Lifting Operations and Lifting Equipment Regulations (LOLER) Examinations. A LOLER examination may only be conducted by a person trained and qualified to do so. Space is provided on the sling label to record LOLER details.



3. DAMAGED OR FAULTY SLINGS



Care & Independence have a strict *Returns & Repair Procedure* to ensure the safety of all involved.

Please do not send slings and bathing products back without completing the Return/Repair form.

4. PRODUCT CARE



4.1 GENERAL CARE:

4.1.1 Repeated washing may lead to sling deterioration. **Care & Independence strongly recommend replacing the sling at the first signs of wear or damage, as opposed to after a specified time frame.**

4.1.2 Each sling is marked with its maximum wash temperature.

4.1.3 Discolouration to some materials or fading labels may occur as a result of laundering.

NOTE: It is a health & safety requirement that the sling label must be legible throughout the sling's working lifetime.

4.2 **WASHING THE SLING:** Every sling is marked with its maximum wash temperature. Close hook & loop fastenings before laundering to prevent damage to the fabric and webbing.

4.3 DECONTAMINATION:

Slings can withstand a strong decontamination routine but regular use may impair the lifespan of the sling.

A copy of these instructions should be provided to the person(s) laundering the sling.



- Wash gently at 40°C. Use a mild soap solution or a non-biological detergent without bleaching additives.
- Do not add bleach or stain removers.
- Air dry, cabinet dry or tumble dry on a cool cycle (below 60°C). Do not tumble dry on a hot setting.
- Do not place on hot pipes or heaters. Keep away from flames.
- Do not iron.

FURTHER INFORMATION: Available to download from our website



- Sling Care Instruction Guide (includes decontamination advice).
- Returns & Repair Procedure



6. FITTING THE SLING

Both COMFORT VALUE™ Recline models are In-Chair slings. The user can be sat on the sling for a period of time following a risk assessment. Also suitable for amputees (single, bi-lateral above-knee). Before you start, check that the sling is in safe and clean condition and is the appropriate size and model for the person to be transferred. There are a variety of ways to use and adjust this sling, so please read carefully.

NOTE: Recline sling models are not suitable for fitting in a chair.

6.1 FROM A SUPINE POSITION

To prevent falls, always roll a person towards yourself, a wall or a bed rail.


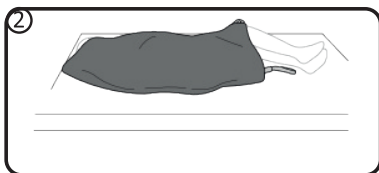
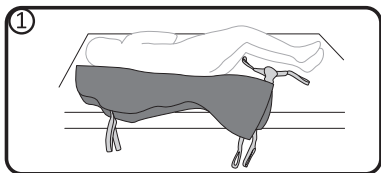


6.1.1

- i. Roll the user onto their side.
- ii. Fold sling in half lengthways with the label on the inside.
- iii. Position the sling with the darts aligned with the user's hip and the centre of the sling in alignment with the user's spine.

6.1.2

- i. Fold the uppermost section over the user like a blanket.
- ii. Tuck the remaining sling material under the user, ensuring straps are tucked in.

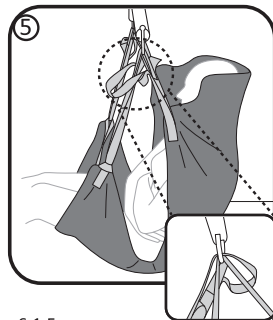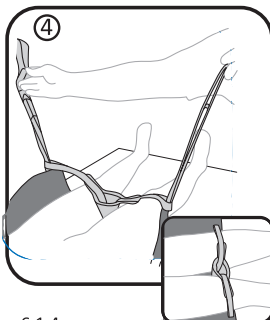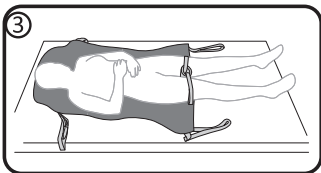

6.1.3

- i. Roll the user onto their opposite side.
- ii. Unfurl the sling beneath the user.
- iii. Reposition user onto their back with sling now spread out fully underneath them.
- iv. Ease the pommel through between the knees.



6.1.4

Thread the leg lift straps through the pommel straps on each corresponding side.

6.1.5

- i. Raise bed to a seated position.
- ii. Attach first leg lift straps then shoulder lift straps to the hoist spreader bar.
- iii. After checking nothing is snagged, complete the transfer.

6.2 HOISTING IN A SEATED POSITION

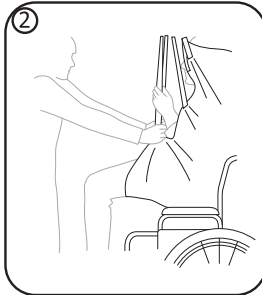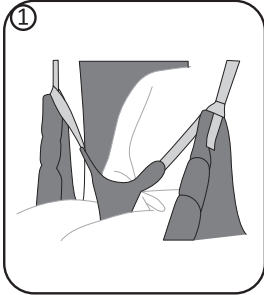

6.2.1 To Hoist **Out** of a Chair:

- i. Thread the leg lift straps through the pommel straps on each corresponding side.
- ii. Attach leg straps to the hoist first then shoulder straps.
- iii. Check for snagging.
- iv. Complete the transfer.

6.2.2 To Hoist **Into** a Chair:

To assist guidance into a chair, a preferred technique is to use your body weight from the front, pushing on the bottom of the leg lift straps, as the user lowers into the chair.

NOTE:

Leg lift and pommel straps can be tucked into corresponding pockets to prevent snagging when not in use.

7. ADDITIONAL SUPPORT AVAILABLE



Free product training is also offered. Please contact us for further information.



- Expert advice
- Full support



SCAN SLING LABEL QR CODE FOR;


- Product information
- User Guide
- Fitting video
- Sling Care Instruction Guide
- Sling Pre-Check Inspection Guide
- Sling Sizing Guide



- Faded label replacement
- Sling repair
- LOLER examinations

5. GETTING TO KNOW YOUR SLING

5.1 THE RECLINE-STYLE SLING

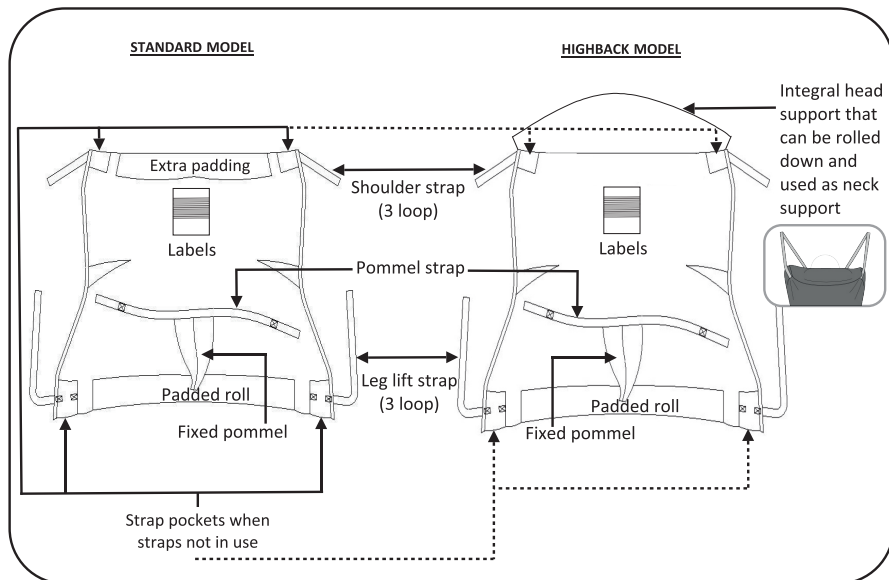



5.2 SLING LABELS EXPLAINED

Each sling has one printed label sewn onto the back.

The label contains useful information such as Sling Model Type, Size, and Serial number. Safety information and washing instructions are included too. There is also a section where LOLER checks can be recorded as well as space to write the User's name.

Further product information can be obtained by scanning the QR code that leads direct to the product webpage where support resources are available such as full product information, user guides, fitting videos and more besides.

Care & Independence QR CODE LEADING TO PRODUCT INFORMATION AND USER GUIDES 		CARE & INDEPENDENCE CONTACT INFO
SLING INFORMATION MODEL SIZE SERIAL NUMBER YEAR OF MANUFACTURE 		LOLER RECORD INPUT
SAFETY INFO		
RECOMMENDED WASHING		SWL IN KG
QUALITY MARKS & INFO		
NAME: _____		

Comfort VALUE™ USER GUIDE

Care & Independence

www.careandindependence.com




Care & Independence Ltd,
Unit 3, Telford Business Centre,
Telford Road, Bicester, OX26 4LD

033 012 80808

sales@careandindependence.com

Monday - Friday 08.30 -17.00

*If a third party supplier has provided the COMFORT
VALUE™ sling, you may like to write their contact
details here:*

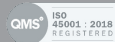

©2023 Care & Independence Ltd All Rights Reserved.
Contents of this publication are intended as a guide only. Although due care and attention has been given to ensure that instructions are accurate and correct, in case of doubt, please refer to your prescribing clinician or contact the Care & Independence team on 033 012 80808. Care & Independence Ltd reserve the right to modify, add or delete details and is under no obligation to notify the user of amendment. All data is correct at time of publication.