

USER INSTRUCTIONS

ECONOMY TOILET FRAME RANGE

| | | | |
|--------|---------------------------------------|--------|-----------------------------------|
| M00869 | Economy Toilet Frame 17" (430 mm) | M87203 | Economy Toilet Frame 19" (480 mm) |
| M00870 | Economy Toilet Frame Width Adjustable | | |

⚠ General Safety Warning ⚠

Please make sure that you read and follow these instructions carefully. Keep these instructions in a safe place for future reference.

Maximum User Weight



1. INTENDED USE

The height and width* adjustable Economy Toilet Frame is designed to assist users when lowering onto and rising from the toilet. It provides comfortable and stable support when getting on and off the toilet and is suitable for use in conjunction with a raised toilet seat. An optional floor fixing kit is also available (M24308).

**M00870 only.*



Figure 1: M00869



Figure 2: M00870

2. SAFETY

WARNING: The maximum user weight must not be more than 190 kg (30 st).

DO NOT: Attempt to stand or climb on the product.

MAKE SURE: All legs are set at an equal height to avoid rocking.

3. UNPACKING

Remove the packaging and check all parts carefully. Contact your supplier immediately if parts are missing or damaged and do not try and assemble.


4. ASSEMBLY

These products are supplied fully assembled.

5. INSTALLATION

Adjust the height and width* for the toilet and user before installation.

5.1 Adjust the Height - All Models

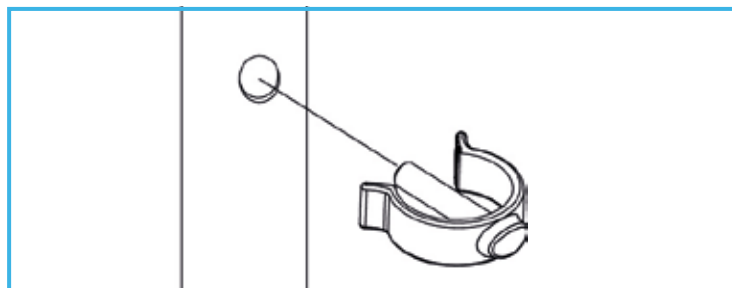

1. Remove the 4 x leg extension E-clips.
2. Slide the leg extensions to the required height, align the holes and replace the E-clips.

Note: Make sure the leg extensions have all been set to the same height.

5.2 Adjust the Width - M00870 only

1. Remove the 1 x cross bar E-clip.
2. Slide the side frame assemblies to the required width, align the holes and replace the E-clip.

Note: Make sure the seat and toilet frame are fully secure before use.

5.3 Fitting the Optional Floor Fixing Kit M24308

Please refer to the user instructions provided with the optional floor fixing kit.

6. MAINTENANCE

Check that all fixing clips are inserted correctly and that there is no excessive movement.

If a floor fixing kit has been fitted, check the fixing screws are secure.

Inspect the product regularly for signs of wear or damage, pay particular attention to the condition of the ferrules.

Any suspected faults or damage should immediately be reported to the supplier and the product taken out of use.

7. CLEANING INSTRUCTIONS

Clean regularly using a general purpose, neutral detergent and warm water.

Please check the suitability of the cleaning products before use. Do not use abrasive cleaners or cloths as these could permanently damage the surface and invalidate any warranty.

To disinfect, wipe with a dilute bleach solution and rinse thoroughly with fresh water.

Allow the product to dry fully prior to use/storage.

8. ACCESSORIES & SPARE PARTS

M24308: Economy Range Floor Fixing Kit

For further information please visit our website or call Customer Services on the details below.

9. PRODUCT DISPOSAL

This product must be thoroughly cleaned/sterilised prior to disposal.

This product is manufactured from metal and plastic and can be disposed of through Local Authority Household Waste and Recycling Centres or through Local Trade Waste Collectors.

10. PRODUCT GUARANTEE

NRS Healthcare guarantee this product for a period of 24 months, from the date of purchase, against faulty materials and workmanship.

This guarantee does not cover abuse, interference or tampering with any part of the products. Issues arising from fair wear & tear and customer misuse are not covered.

This guarantee covers domestic use only. If the product is used in commercial or semi-commercial settings such as care homes, schools, rental or communal accommodation etc, the guarantee is valid for one year only. Consumable parts are also guaranteed for one year only.

Any guarantee issued by NRS Healthcare does not affect your statutory rights.

NRS HEALTHCARE Sherwood House, Cartwright Way, Forest Business Park
Bardon Hill, Coalville, Leicestershire, LE67 1UB, UK.

Tel: Customer Services: +44 (0)345 121 8111

Email: customerservice@nrshealthcare.co.uk

NRS Healthcare is a leading supplier of community equipment and daily living aids.
Please visit one of our websites to view our full range of products:
www.healthcarepro.co.uk or www.completecareshop.co.uk