



# Apolo

### Apollo 5 Premium

Product Code - APH183



Max Weight Limit: 28 stone / 177.8 Kg Min Weight Limit: 6 stone / 40 Kg **Risk Factor:** Medium to High Risk Warranty: 2 years Cell Design: 18 individual cells **Mattress Dimensions (inflated):** 200cm x 88cm x 12.5cm Mattress Weight: 3.5 Kg **Power Unit Dimensions:** 26cm x 12cm x 10cm Power Unit Weight: 2 Kg

#### Electrical

Safety Standards: BS EN 60601-1 BS5724-1 (CE, FDA, TUV) Voltage: AC 220v-240v 50Hz/60Hz Noise Level: below NC27 Power Input: 12w Fuse: T1A, 250v



Apollo Healthcare Technologies Limited, Holme Street, Liversedge, West Yorkshire WF15 6JF Tel: +44 (0) 1924 614567 Fax: +44 (0) 1924 607480 Email: sales@apollo-ht.co.uk www.apollo-ht.co.uk

## **Dynamic Overlay Mattress**




#### www.apollo-ht.co.uk

### Apollo 5 Premium<sup>™</sup> Dynamic Overlay Mattress



The **Apollo 5 Premium** dynamic overlay system is suitable for those considered at '**Medium to High Risk'** of developing pressure ulcers, offering optimum pressure relieving qualities and comfort. The simple to use pump ensures this ever popular overlay can be used in most nursing or homecare environments.

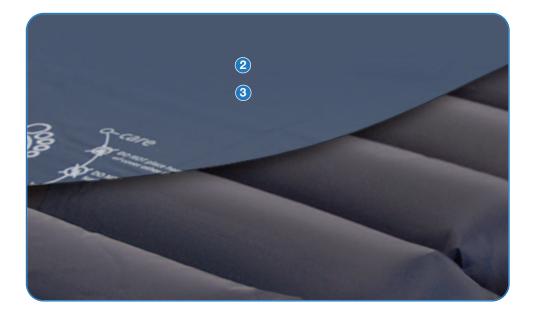

- 28 stones weight limit
- Medium to High Risk
- Overlay system
- Manual adjustment to user's body shape and position for optimum pressure care
- Designed for users at medium to high risk, or have existing tissue damage (pressure ulcer grade 1 to 4)
- Dynamic therapy alternating action for lower interface pressures and gentle capillary stimulation
- Humidity management to help prevent tissue breakdown and skin maceration
- Enhanced shear and friction protection
- Zipped cover is easy to remove for easier cleaning
- Clinically effective against bacteria/viruses, with proven resistance to MRSA
- Fire retardant to Standards BS7175 Crib 5
- Clear and simple controls reduces need for complex training
- Removable cells for lower cost servicing and repair
- The mattress cover can be wiped down or completely removed and laundered up to 95°C to aid decontamination protocols
- 1 Adjustable weight settings
- 2 PU cover to prevent fluid ingress and internal mattress contamination
- 3 Waterproof with high vapour permeability for improved infection control
- 4 Quick access CPR for emergency mattress deflation

